

For Veterinary Use Only

READ ALL INSTRUCTIONS BEFORE BEGINNING THE TEST

RIDX™ ASFV Ag Test Kit

[Catalogue Number: LGM-PAG-11]

Introduction

African swine fever virus (ASFV), a large double stranded DNA virus, genus *Asfivirus* in the only member of the *Asfarviridae* family¹, is known as not infectious to humans and not directly affect public health².

African swine fever (ASF) is a lethal hemorrhagic fever of pigs with high mortality and morbidity rates approaching 100 percent³. ASF was first reported in Southeast African countries, including Kenya, in the 1920s³. The spread of ASF outside Africa where was first described is persistent in Europe, the Russian Federation, China, Mongolia and Vietnam. It also has posed a serious threat of animal health, pig production and global food security. Since its introduction to and dramatic spread in 2018 in China, which has half of the world's pig population, current epidemiological situation keeps impact on the pig industry⁵.

Apparent symptoms of disease are eruptions, in which the skin turns pinkish and nearly purple marks on the skin of the extremities, ears, chest, abdomen and perineum. Moreover, abortion can frequently occur in gestating females and be sometimes the first clinical sign of an outbreak².

At present, there are no treatment or vaccines available to prevent ASF so the control strategy mainly relies on strict administration of sanitary measures and slaughtering of infected and exposed animals. Therefore, rapid and accurate diagnosis of the infections is essential for prevention, control, and extermination of the disease in affected various countries⁶. Point-of-care rapid and sensitive detection of ASFV is key to the timely implementation of control⁷.

Principle

The RIDX™ ASFV Ag Test Kit is a lateral flow chromatographic immunoassay for the qualitative detection of ASFV in porcine blood.

This kit shows two letters which are the test (T) line and the control (C) line on the surface of the device. If the ASFV antigen exists in the sample, it binds to the gold-conjugated anti-ASFV antibody. The antigen-antibody complex moves through the membrane by capillary force and responds to the secondary anti-ASFV antibody on the test line, resulting in a red line. The control line indicates that the test is performed correctly and should appear when the test is complete.

Two different monoclonal antibodies to the conservative phosphoprotein p30 protein⁸ of ASFV are used as a capture and detector in the kit. The RIDX™ ASFV Ag Test Kit can detect ASFV in porcine blood with high accuracy.

Performance

1. Sensitivity & Specificity

		PCR		Total
		+	-	
RIDX™	+	10	0	10
ASFV Ag	-	2	435	437
Test	Total	12	435	447

Sensitivity: 83.33% (10/12, *95% CI: 55.20% ~ 95.30%)

Specificity: 100% (435/435, 95% CI: 99.12% ~ 100%)

Diagnostic Agreement: 99.55% (445/447, 95% CI: 98.38% ~ 99.88%)

* 95% CI: 95% Confidence Interval

2. Limit of Detection: 5×10^4 TCID₅₀/mL, 7.5 HAD₅₀/mL

3. Cross-Reactivity

Potentially cross-reactive substances listed below have no effect on the performance of the RIDX™ ASFV Ag Test Kit.

Pathogen	Titer (TCID ₅₀ /mL)
Classical swine fever virus (CSFV)	1×10^5
Japanese encephalitis virus (JEV)	1×10^5
Porcine circovirus type 2 (PCV2)	1×10^4
Porcine reproductive and respiratory syndrome virus (PRRSV)	1×10^3
Porcine epidemic diarrhea virus (PEDV)	1×10^4
Transmissible gastroenteritis virus (TGEV)	1×10^3

Kit Components


Component	Number/Kit
1 ASFV Ag Test device	10
2 Anticoagulant tube with EDTA	10
3 Disposable capillary tube	10
4 Dilution buffer (4 mL)	1
5 Instructions for use	1

Storage & Stability

1. Store the test kit at 2~30°C (35.6~86.0°F). **Do not freeze.**
2. Do not store the test kit in direct sunlight.
3. The test kit is stable within the expiration date marked on the package label.

Sample Preparation


1. Use **swine whole blood** as a specimen.
2. Collect 1 mL (0.5~1.5 mL) of the whole blood sample and put it into an anticoagulant tube immediately.
3. Close the cap on the anticoagulant tube and invert the tube 5 times to mix blood sample and ethylene diamine tetra acetic acid (EDTA).

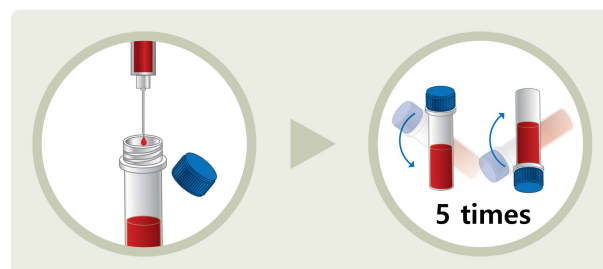

4. The anticoagulated whole blood samples should be used instantly after collection. If you cannot use the samples immediately, store them refrigerated (2~8°C/35.6~46.4°F) or keep them on ice. **Do not freeze the anticoagulated whole blood samples.**
5. For the best performance, the blood samples should be tested within 24 hours after collection. But, results from stored samples may differ from those obtained before storing.
6. When using stored anticoagulated blood, use only blood that has been stored for up to 3 days (72 hours) at 2~8°C (35.6~46.4°F), a refrigerated condition. **Do not use frozen samples.**

◆ Test Procedure

1. All samples and test components should be at room temperature (15~30°C/59~86°F) before use.
2. Take 20 µL of the anticoagulated whole blood to the black line of the capillary tube.
3. Apply the taken sample into the sample hole (S).
4. After 1 to 3 minutes, which is the time for the blood sample to be sufficiently absorbed into the sample pad, apply 3 drops (~100 µL) of the dilution buffer into the sample hole (S).
5. Read test result at 15 minutes. **Do not read results that appear after 15 minutes.**

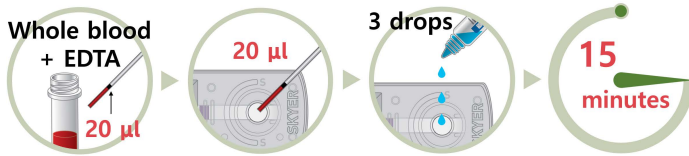


[Summary of Test Procedure]

◆ Interpretation of Results

1. Positive result


Test (T) line and control (C) line within the result window indicate the presence of African swine fever virus antigens.

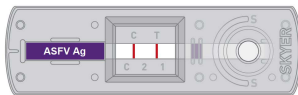

※ If the blood specimen is too viscous or hemolyzed, flow along the membrane may be impeded, resulting in nonspecific false-positive results. Therefore, the results observed after the designated time are deemed unreliable.


2. Negative result

Only control (C) line appears in the result window.

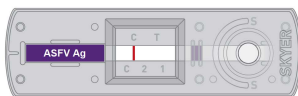

3. Invalid results

If the control (C) line does not appear, the result might be considered invalid. The sample should be retested.

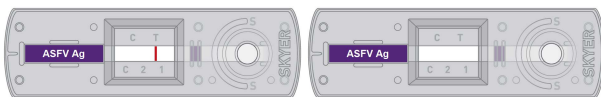

◆ Precautions

1. This test kit is for veterinary *in vitro* diagnostic use only for pigs. Do not use this test kit for other animals.
2. This rapid kit is only for preliminary screening. The final decision should be made by a qualified veterinarian based on the results of this kit, clinical symptoms and evaluation by a veterinarian, and, if necessary, the results of additional detailed diagnostic procedures.
3. The test device is sensitive to humidity and heat. Use the test device within 10 minutes after removing the foil pouch.
4. Do not touch the membrane in the sample hole on the device.
5. The device should not be used if the foil pouch is damaged.
6. Do not use an expired test kit. The expiration date is marked on the package label.
7. Do not reuse the test components.
8. Do not mix components from different lot numbers because the components in this kit have been quality control tested as a standard batch unit.
9. Decontaminate and dispose of all samples, used kits, and potentially contaminated materials following national and local regulations.
10. All samples should be handled as being potentially infectious. Wear protective gloves while handling specimen and used components. Wash hands thoroughly afterward.

◆ References

1. International Committee on Taxonomy of Viruses (ICTV). *Virus Taxonomy: 2020 Release*. Email Ratification March 2021 (MSL # 36).
2. Sanchez-Vizcaino JM, Laddomada A, Arias M. African Swine Fever Virus. *Diseases of Swine*, Eleventh Edition, Eds: Zimmerman JJ, Karriker LA, Ramirez A, Schwartz KJ, Stevenson GW, Zhang J. 2019; 25: 443–451. John Wiley & Sons, Inc., NJ, USA.
3. Costard S, Wieland B, de Glanville W, Jori F, Rowlands R, Vosloo W, Roger F, Pfeiffer DU, Dixon LK. African swine fever: how can global spread be prevented? *Phil Trans R Soc B*. 2009; 364: 2683–2696.
4. Penrith ML, Thomson GR, Bastos ADS. African swine fever. *Infectious diseases of livestock* (Eds: Coetzer JAW, Tustin RC). 2004; vol 2: 1088–1119. Oxford University Press, Oxford, UK.
5. Dixon LK, Sun H, Roberts H. African swine fever. *Antiviral Res*. 2019; 165: 34–41.
6. Sastre P, Gallardo C, Monedero A, Ruiz T, Arias M, Sanz A, Rueda P. Development of a novel lateral flow assay for detection of African swine fever in blood. *BMC Vet Res*. 2016; 12: 206.
7. Wang X, Ji P, Fan H, Dang L, Wan W, Liu S, Li Y, Yu W, Li X, Ma X, Ma X, Zhao Q, Huang X, Liao M. CRISPR/Cas12a technology combined with immunochromatographic strips for portable detection of African swine fever virus. *Comm Biol*. 2020; 3: 62.
8. Afonso CL, Alcaraz C, Brun A, Sussman MD, Onisk DV, Escibano JM, Rock DL. Characterization of P30, a Highly Antigenic Membrane and Secreted Protein of African Swine Fever Virus. *Virology*. 1992; 189: 368–373.

◆ Symbol Descriptions

	License number
	Catalogue number
	Batch code, Lot number
	Consult instructions for use
	Contains sufficient for <n> tests
	Do not reuse
	<i>In vitro</i> diagnostic medical device
	Temperature limitation
	Do not use, if the package is damaged
	Upper side
	Manufacturer



SKYER, INC.

#532, 416, Hwagok-ro, Gangseo-gu, Seoul, 07548,
Republic of Korea
TEL: +82-2-706-6801, FAX: +82-50-4096-6988
Technical support: marketing@skyer.co.kr
www.skyerdiagnostics.com

Korean Veterinary Diagnostics Manufacturer License No. 300